



Installation Instructions (Linux/Windows)-3DSlicer EPCMR module (preRelease 08-05-2026):

1. Download and install Slicer (any version 5.7+ will do)
2. Download the python module containing .zip file from polybox:
<https://polybox.ethz.ch/index.php/f/4515989317>
3. Unzip EPCMR module to a location of your choice; at the root level there is the main module file named "EPCMR.py". In your file explorer, drag and drop this EPCMR.py file on 3DSlicer application window.

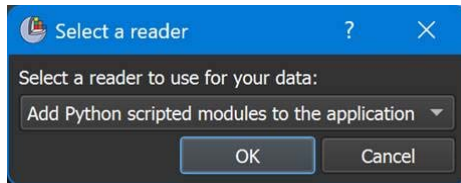

"Select a reader to use for your data: - Add Python scripted modules to the application ?" - press OK - confirm in the next popup window "Add selected module to additional module paths ?" - "YES" ==> **DONE !**

4. Restart 3DSlicer app - our EPCMR logo should appear in the upper module toolbar - press it and enjoy !
5. You may want to set your Slicer application to use dark mode settings (one-time config):

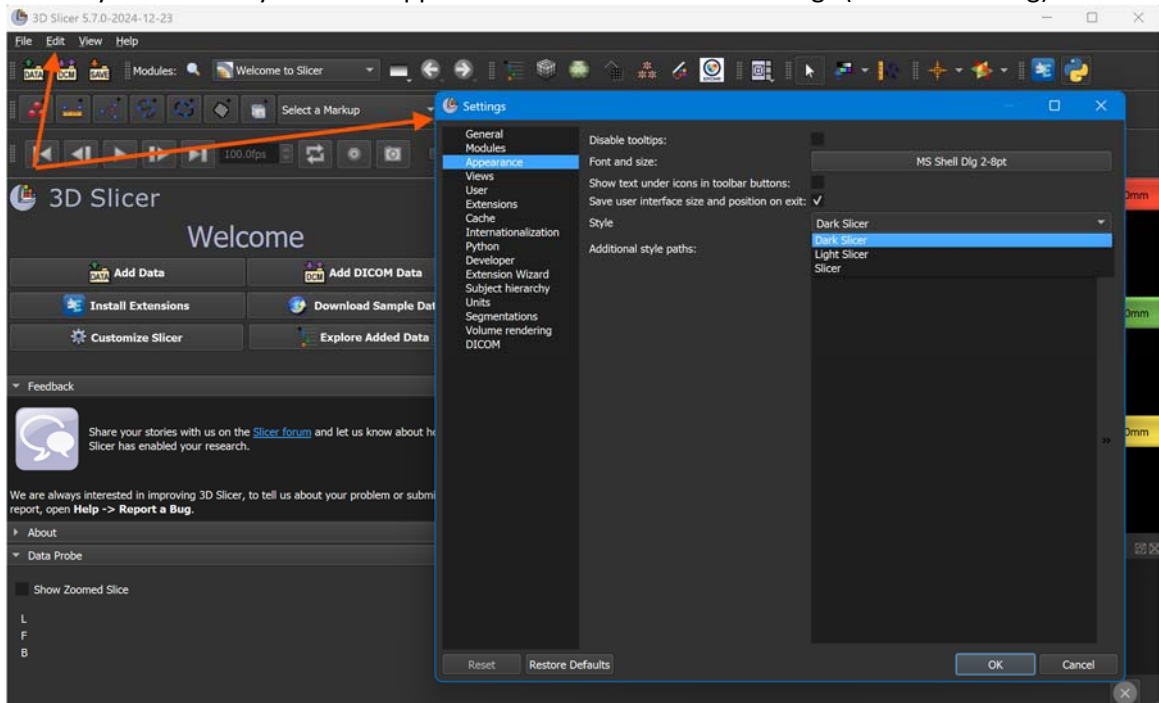

6. Module root level of EPCMR folder contains "**ExampleData**" folder, use as follows:
 - a. **stlModels**: simply drag and drop the .STL files into the slicer app window (models are loaded/colored etc. immediately as per our implemented visual style)
 - b. **readPickleShotList**: simulate OpenIGTLink live transmission of transformMessages (= "CMR scanner runs and delivers transform messages...")
 - i. simply run in your python interpreter "ReadPickle_shotList09.py" - it will read the data from same level pickle file named "tracking_data_2025-11-06_13-35-02_repaired.pkl"
 - c. **replay_jsonlines**: in the running EPCMR module you can open the "igt_l_stream.jsonl file" under "CatheterReplayer" (colored upper bar can be clicked to toggle "LIVE" and "REPLAY" mode; only during "REPLAY" you can pick "Open .jsonl file" and load it...

IMPORTANT NOTE:

Your standard 3D Slicer download does *not* include the required OpenIGTLink Python library by default. Hence, on its first run, the EPCMR module automatically installs the pyigtl==0.3.1 dependency from this source: <https://pypi.org/project/pyigtl/>